

PCT / CA 98 / 00249

23

APRIL 1998

(23-04-98)

REC'D 03 JUN 1998

WIPO

PCT

THE UNITED STATES OF AMERICA

TO ALL TO WHOM THESE PRESENTS SHALL COME:

UNITED STATES DEPARTMENT OF COMMERCE
United States Patent and Trademark Office

April 7, 1998

THIS IS TO CERTIFY THAT ANNEXED HERETO IS A TRUE COPY FROM THE RECORDS OF THE UNITED STATES PATENT AND TRADEMARK OFFICE OF THOSE PAPERS OF THE BELOW IDENTIFIED PATENT APPLICATION THAT MET THE REQUIREMENTS TO BE GRANTED A FILING DATE UNDER 35 USC 111.

APPLICATION NUMBER: 60/041,280

FILING DATE: March 21, 1997

By Authority of the
COMMISSIONER OF PATENTS AND TRADEMARKS



L. BLAND

Certifying Officer

69/041280

PROVISIONAL APPLICATION COVER SHEET

Request for filing a PROVISIONAL APPLICATION under 37 CFR 1.53(b)(2).



DOCKET NUMBER		1960.15 PV		Type a plus sign (+) inside this box.	+
INVENTOR(s)/APPLICANT(s)					
LAST NAME	FIRST NAME	MIDDLE INITIAL	RESIDENCE (CITY AND EITHER STATE OR FOREIGN COUNTRY)V		
1. RICHARDSON	1. TIMOTHY	M.	1. 16 MATSON DRIVE BOLTON, ONTARIO, CANADA		
2.	2.	1.	2.		
3.	3.	2.	3.		
TITLE OF THE INVENTION (280 characters max)					
MICROSCOPE SLIDE SYSTEM AND METHOD OF USE THEREOF					
CORRESPONDENCE ADDRESS					
FITZPATRICK, CELLA, HARPER & SCINTO 277 Park Avenue New York, N.Y. 10172 Telephone No. (212) 758-2400					
ENCLOSED APPLICATION PARTS (check all that apply)					
<input checked="" type="checkbox"/> Specification	Number of Pages	21	<input type="checkbox"/> Small Entry Statement		
<input checked="" type="checkbox"/> Drawing(s)	Number of Sheets	12	<input type="checkbox"/> Other (specify)		
METHOD OF PAYMENT (check one)					
<input type="checkbox"/> A check or money order is enclosed to cover the Provisional filing fees			PROVISIONAL FILING FEE AMOUNT (\$)		\$150.00
<input checked="" type="checkbox"/> The Commissioner is hereby authorized to charge filing fees and credit Deposit Account Number: 06-1205					

The invention was made by an agency of the United States Government or under a contract with an agency of the United States Government.

☒ No.☐ Yes, the name of the U.S. Government agency and the Government contract number are: _____

Respectfully submitted,

SIGNATURE

Richard P. Bauer

DATE March 21, 1997

TYPED or PRINTED NAME RICHARD P. BAUER

REGISTRATION NO. 31,588
(if appropriate)

☐ Additional inventors are being named on separately numbered sheets attached hereto.
PROVISIONAL APPLICATION FILING ONLY

F501\W133317\dcw

50043588-03437



60/041280 ✓

MICROSCOPE SLIDE SYSTEM AND METHOD OF USE THEREOF

3

4 Field of the Invention

5

6 The present invention relates to a microscope slide system and a method of
7 preparing specimens on microscope slides.

8

9 Description of the Prior Art

10

11 The practise of making well-prepared microscope slides has been accomplished
12 in the past by mounting samples in a mounting media and fixing a cover slip to the
13 mounting media by the adhesive action of the mounting media or with sealing rings of
14 another adhesive material. The process of cleaning and sterilizing the slides,
15 maintaining this sterile condition while handling them during the mounting process and
16 applying the mountant and the sample involved a great deal of skill and attention. This
17 process was even more demanding when slides were prepared for use in very high
18 magnification, or polarized light or in dark field illumination since tiny contaminating
19 particles, irregularities or scratches in the slide, mountant, sample or cover glass could
20 obscure critical parts of the sample image or could be misinterpreted as either part of
21 the sample or conversely a real part of the sample could be misinterpreted as an
22 "artifact".

23

24 The preparation of live samples for real time living microscopic examination or
25 the culture of cells, bacteria, viruses or other living materials presented technical
26 problems in the preparation of the slide, both from a sealing perspective and from a
27 hazard protection standpoint.

28

/61220 0214003

1 It is an object of the present invention to provide a novel microscope slide
2 system. It is a further object of the present invention to provide a novel method of
3 preparing specimens on microscope slides

4
5 **Background of the Invention**
6

7 This present invention is intended to address the needs of the routine and
8 research microscopist by providing a simple and effective slide system and method of
9 producing high quality slides which are free of contamination and "artifacts" and
10 maintain sterile conditions. These slides also provide a range of features which will aid
11 the researcher in using the slides to obtain data about the sample by supplying reference
12 markings, filtering the light reaching the sample, and providing means to vary the
13 environmental conditions in the sample space. The slides are designed to provide
14 optimal conditions for live samples and for the long term maintenance of life supporting
15 conditions in the sample space. Further features provide safeguards to reduce the risk
16 of hazardous materials in the sample space moving out of the sample space and reaching
17 the working environment.

18
19 **Detailed Description Of The Invention**
20

21 In its simplest embodiment, the present invention provides a disposable, single
22 use, sterile slide system which offers quick preparation, use and disposal of the slide.
23 In its most basic form the slide system comprises a standard microscope slide but with
24 an adhesive layer applied to the slide. A microscopy sample is placed within the area
25 contained by the adhesive layer. The cover slip is then applied to the adhesive layer to
26 form a sealed space for the sample. The adhesive layer secures the cover slip to the
27 slide, preventing the sample from escaping, or the solvent(s) in the sample from
28 evaporating or being contaminated and preventing the cover slip from moving under the
29 influence of immersion oil. The system of slide and cover glass forms a stable and easy
30 to handle unit once sealed. The final thickness of the sample space is established by the

202 317 8136-7 6/96

1 adhesive layer thickness, the geometry of the base slide and the closing force applied to
2 the cover slip during the cover slip application and adhesive activation operation. In
3 other embodiments, a seal member such as a neoprene, Delrin or glass ring element can
4 be positioned between the cover slip and the slide to establish as desired final thickness,
5 with adhesive maintaining the seal member between the cover slip and the slide.
6 Integral expansion areas in the sample space allow trapped gas to act as an absorbing
7 space for excess sample material during the closing operation.

8
9 The adhesive layer can be applied to either the slide, the cover slip or both,
10 depending on the application. For instance in cases where cells are cultured on the
11 cover glass, the adhesive layer would only be coated on the slide. The particular
12 adhesive employed is not particularly limited and any suitable adhesive as will occur to
13 those of skill in the art can be employed. For example XYZ as manufactured by ABC
14 can be employed.

15
16 The adhesive can be applied as a one part system coated onto the slide or it
17 could be applied as a precoated and diecut double sided tape form. The adhesive could
18 also be supplied as a two part system with a first part applied to the slide and an
19 activating second part applied to the cover glass, in operation the adhesive bond would
20 only be formed by the catalytic or stoichiometric reaction of the two adhesive parts.
21 The adhesive could be a permanent hardening adhesive such as 3M automotive trim
22 adhesives, or it could be a long term flexible and / or removable adhesive such as those
23 used in adhesive bandages or tapes.

24
25 Where the ability to remove the cover slip and rework the sample, or add to it,
26 is desired, a low tack adhesive could be used on the cover slips or slides. These cover
27 slips could be changed at will to add or remove or work with materials in the sample
28 space. In the case of the some of the adhesives mentioned above, once coated, the
29 adhesive would be protected by a removable peel off protecting material or film. In
30 order to prepare the adhesive for use the film would be peeled away to prepare the

150043280-032409

1 adhesive to contact the desired surface. The adhesive layer thickness can be tightly
2 controlled during coating to ensure that the final closed sample space thickness is
3 correct for the particular slide application. Different adhesive thickness can be used to
4 provide a controlled thickness spacer to ensure high accuracy and different sample
5 thicknesses.

6
7 The adhesive can be controlled for refractive index and scattering properties to
8 allow the adhesive area to become a source of side lighting for dark field or
9 ultramicroscopic applications. The adhesive layer can contain phosphors or fluorescent
10 compounds to act as a side light source. The phosphors and fluorescent compounds
11 placed in the adhesive compound can also act as quality control features during
12 application, use and storage. In this way, the adhesive can be checked under UV light
13 for gaps, voids or contamination after coating or after cover slip sealing.

14
15
16 The slide and cover can be made of any material that is transparent at the
17 wavelengths of interest. For visible light use, optical plastics, glasses such as BK7 or
18 other suitable materials can be used, as will occur to those of skill in the art. For IR
19 use materials such as silicon or AMTIR can be used or other suitable materials. For
20 ultraviolet light materials such as fused quartz, crystalline quartz, sapphire, spinel,
21 zircon, diamond, calcium fluoride, lithium fluoride, or magnesium fluoride can be used
22 or other suitable materials. The cover material can be chosen from materials such as
23 those listing above and can be the same material as the base slide or can be chosen to
24 match, complement, or correct for, the refractive index and dispersion of the sample or
25 mounting media or for the immersion fluid. For deep UV and X-ray use, thin metallic
26 films and crystalline films can also be used for slide and cover slip.

27
28 The base slide can be in conventional 1 inch by 3 inch form or can be round,
29 square or geometrically patterned to suit the application.

30

60041280.032197

1 The slide material in the sample area can be etched to give diffuse lighting of the
2 sample. The optical space may also be tapered by etching, machining or any other
3 suitable technique to give various depths of sample and to assess particle or entity size
4 of the sample or sample components. The base slide can also include polarizing
5 material to give polarized light in the sample space.

6
7 The material forming the slide may be coated with interference layers or may be
8 ionically coloured to produce monochromatic or polymonochromatic illumination at the
9 sample space via a novel process. This principle can be used to form a low cost three
10 colour translation system used with polychromatic illumination and three stage post
11 filtering to produce a visible translation of an ultraviolet, visible and / or infrared
12 image. The process of colour translation involves the use of three or more
13 monochromatic light sources to create three monochromatic images which are then
14 combines to create a full colour image. In this version of colour translation where the
15 slide provides a filter system incorporated into the slide surface using thin film
16 interference filters, the light from a broadband source is filtered into three
17 monochromatic bands by the slide and the final separation of the three bands into the
18 three separate monochromatic images carried out using filters between the microscope
19 objective and the imaging system. This process allows the microscope to be used as a
20 colour translating system without needing a monochromator in the light path before the
21 slide. A typical use would be to convert 256, 273 and 282 nanometer wavelength light
22 into a corresponding blue, green and red television images.

23
24 The base slide can also contain more than one layer of transparent material to
25 correct for aberration or allow conduction of light along the XY plane of the slide by
26 total internal reflection.

27
28 In a normal microscope, focusing of the objective on the slide is a matter which
29 requires careful adjustment of the fine focus and coarse focus of the microscope and
30 may involve some degree of difficulty for people who are not experienced with

E0041230-032197

1 microscope operation. It can often also be the source of damage either to the objective
2 or to the slide caused by the objective crashing into the slide surface and damaging the
3 objective or the slide surface.

4

5 It is an advantage of the present invention to allow the possibility of this
6 happening to be reduced and, if the slides are prepared properly using the system
7 described in this patent, to substantially reduce the possibility of the objective crashing
8 into the slide by carefully controlling the closed thickness of the slide cover adhesive
9 and slide system. The method of controlling the closing thickness is to use specific
10 thicknesses of adhesives and/or spacers and tightly controlled thicknesses of cover
11 glasses and tightly controlled thicknesses of slides.

12

13 When the resulting slides are placed in a holder of carefully controlled
14 thickness, the resulting z-axis location of the sample space and the upper surface of the
15 cover slip make it possible to quickly move the objective to this known reference
16 position under computer driven control with no risk of crashing the objective into the
17 cover slip. This allows the microscopist to quickly find the proper focus in the center
18 of the sample space and to adjust the focus from this point higher or lower into the
19 sample.

20

21 In some circumstances, it may not be possible or commercially feasible to obtain
22 a desired thickness of adhesive to accommodate the thickness of sample desired for the
23 microscopical application in question. For example, if a sample thickness of 10
24 microns is desired, it might be difficult to find an adhesive which could be applied
25 effectively in a thickness of 10 microns. In such cases, it is possible to chemically,
26 mechanically or otherwise etch back or grind back the slide in the area where the
27 adhesive will be applied so that a thicker adhesive layer can be applied without
28 compromising the final closed thickness of the system. The process of etching back or
29 grinding back a ring for the adhesive can also serve to assist with alignment and
30 location of the slide during the closing process.

60041280-03214009

1
2 An advantage of another embodiment of the present invention is the provision of
3 one or more expansion zones or expansion rings or expansion wells in the sample area.
4 As the cover slip is closed over a small sample with liquid components placed in the
5 XY plane and in the center of the sample area, the liquid will flow outward toward the
6 adhesive zone of the slide. If too much sample material was placed on the slide and
7 there is no expansion area for this material, it will breach the adhesive or the cover slip
8 during the closing operation due to the incompressible nature of the liquid. If surge and
9 expansion areas are provided, where excess liquid can be trapped and air or other
10 ambient gas can be compressed as the liquid flows during the closing operation, then
11 pressure build up can be accommodated or controlled and adhesive seals and cover glass
12 will remain intact. The fluids and gases in the sample chamber of the slide may expand
13 due to heat, illuminating light or chemical reaction after the seal is completed during
14 slide preparation, which in conventional slides can cause the cover slip or the seal to be
15 breached. The expansion zones also provide an element of defense against this form of
16 expansion in use.

17
18 The expansion zones can be in the form of rings or moat in the surface of the
19 slide, wells in the surface of the slide, tapered zones or any other suitable structure as
20 will occur to those of skill in the art. In general, the expansion zones would be outside
21 optically used areas of the sample area. The expansion zones may be created by
22 chemical or mechanical etching of the slide material or by any other suitable means.
23 Representative examples of such means can include ion etching, waterjet or solid
24 particle jet erosion, casting, pressing or grinding of the slide material.

25
26 The expansion zones can either be on the outside surface of the slide or may be
27 formed by burying an expansion zone in a center layer of a two layer slide with
28 connecting ports to the sample space. Also the expansion zones can be in the form of
29 graded wedge shapes in the surface of the slide in the sample area contained in the
30 adhesive ring, to allow for different thickness of sample in the closed area or to allow

E0044380-033197

1 cells to move in different ways depending on the relative thickness or to allow
2 spectroscopic measurements to be made at different sample thickness.

3 Because of the unique hazards presented by working with live cells or live
4 cultures of bacteria or viruses or other biohazardous materials, it's important to make
5 sure that the material both cannot escape from the sample area, and this is accomplished
6 by the adhesive rings, but it is also important to make sure that if material tried to
7 escape from the area it is effectively neutralized. In the present invention,
8 neutralization or the rendering safe of the material is accomplished by coated rings,
9 either in addition to the adhesive rings or as a component within the adhesive rings, of
10 antibiotic, antiseptic or other neutralizing material which for instance could take the
11 form of a radioactive ring.

12
13 These antibiotic or chemical rings can be deposited onto the slide either inside
14 the main adhesive ring or between two concentric adhesive rings or as part of the
15 adhesive rings. In the case wherein the antibiotic or antiseptic is used in the slide
16 design, it can be important to isolate the antibiotic or antiseptic chemical from the
17 sample material itself and this can be accomplished simply by the proper placement of
18 the chemical or antibiotic ring outside the diameter of the expansion ring or mote.
19 Thus the expansion ring serves as a isolation area where a physical boundary is
20 maintained between the antiscptic and antibiotic ring and the sample itself. The
21 problem here would potentially be that if the sample material contacted the antibiotic or
22 antiseptic it might be adversely affected, giving false results of the normal life cycle of
23 the living system being viewed.

24
25 Because the of unique nature of these slides it may be important to present data,
26 scales or shapes right in the sample space for purposes of comparison or record keeping
27 or scaling. It is possible to using chemical methods or laser methods or ion etched
28 methods or ion deposition to create within the sample space reticules or serial numbers
29 or scale bars or reference shapes by which the sample can be compared to known sizes
30 and geometries or where the serial numbers of the sample can be accurately tracked. It

50041380 032197

1 is also possible to record the time that the sample was prepared using a laser etch to
2 time or date stamp or code the slide just prior to the slide closing process, and this
3 implementation of this aspect of the provision of these marking just prior to slide
4 closing will be a subject of another patent related to the slide closing system.

5
6 The features or wells discussed above, in terms of expansion wells, can be
7 created during this same step using the chemical laser or ion etching to etch the wells at
8 the same time as creating the reticules, scale bars or serial numbers.

9
10 As it is often important to test the effect of different substances on samples in a
11 living state and this can be easily accomplished using the slides in accordance with the
12 present invention. This can be accomplished by using known printing or thin film
13 coating techniques to apply thin films of test chemicals, biological agents, antibiotics or
14 antibacterial agents. It can also be accomplished by forming wells in the slide surface
15 and filling or otherwise applying the desired materials to these wells and/or by any
16 other suitable means as will occur to those of skill in the art.

17
18 These methods can also be used to provide the coating of stains, fluorochromes,
19 or immunologically based fluorochromes, commonly known as immune antigen
20 fluorochromes.

21
22 If for instance, a set of six common antibiotics who were coated into the test
23 space and a living sample of a bacteria was put in the test space, it is a simple matter to
24 study the susceptibility or the bacterial to the different test antibiotics by looking
25 through the microscope and seeing which areas of the living sample were adversely
26 affected by the antibiotics which were deposited on the slide. In this way the sterile
27 pre-packaged slides can be opened and have a sample of bacteria placed in the sample
28 chamber with an appropriate nutrient, the whole system will be sealed and then placed
29 under the microscope to watch how the bacteria developed.

30

60041280, 032197

1 This patent also presents a rapid method of testing synergistic effects between
2 different chemicals and antibiotics. For instance in the case of H-pyloria where it was
3 found that a bacteria was not susceptible to standard antibiotics such as tetracycline
4 unless it was also in the presence of a bismuth compound, tests of this type can be
5 rapidly conducted using this embodiment of the present invention. In such cases, an
6 antibiotic can be coated alongside another compound or in conjunction with another
7 compound or mixed with another chemical compound or chemical antiseptic or second
8 antibiotic. The exact effects of the combination of materials or chemicals on the
9 bacterial system can then be studied easily. Using this kind of test gives a very rapid
10 way of doing assays of bacterial characteristics wherein it is not required to grow the
11 bacteria in culture, which may take days or at least one day or more. This gives the
12 clinician a faster way of assessing the biological hazard or drug resistances and
13 susceptibility of a specific bacteria or infectious agent. These methods can be used
14 with viruses, cells, parasites and infectious agents other than bacteria.

15
16 In the case of stains or fluorochromes, vital stains or vital fluorochromes, can be
17 coated or otherwise located on the slide and once the living sample in a suitable media
18 is spaced in the sample space, over the period of a few moments after the sample
19 contacts the stains or fluorochromes, the fluorochromes will dissolve into the media and
20 will be taken up selectively by the applicable portion of the sample. This means that the
21 clinician does not have to do the separate staining step that is normally associated with
22 the vital stained or vital fluorochrome slide.

23
24 Using traditional semiconductor manufacturing techniques it is possible to
25 deposit grids or conductors within the sample space that extend out beyond the adhesive
26 ring or underneath the adhesive ring out to connector pads on the edge of the slide.
27 These grids can be used for electrical testing of living samples or mobile samples within
28 the sample space. For instance, electrophoresis reactions can be carried on under the
29 view of the microscope within the sample space either in reflected light or transmitted
30 light. This allows the researcher to examine precisely electrochemical or bioelectrical

60044280-032197

11

1 effects on living systems or on mobile material samples. This also opens new
2 possibilities in the area of electro-chemistry for studying the precise reactions that are
3 involved in electrochemical deposition or decomposition and interactions.

4
5 Small circuit mount inductors or other components in the slide body can also be
6 placed in the sample area via known surface circuit mount techniques. These can then
7 be coated with suitable adhesives or optical compounds or optical fluids which will keep
8 the sample space intact while still allowing the SMT inductors and their related
9 connection grids to function and set up magnetic fields or alter environmental
10 conditions such as temperature within the sample area. In much the same way as the
11 electrical system mentioned above, these inductors can create fields within the sample
12 space to study the micro effect of magnetic fields on living systems.

13
14 Inductors can also be created by depositing helical, or grid or linear patterns of
15 conductor material, again using semiconductor techniques and to deposit these patterns
16 strategically on the slide in the sample area to create magnetic or electrical fields as
17 desired. These grids can be overcoated with thin layers of inert optically transparent
18 material or with insulating material or with dielectric material or with semiconducting
19 material. All these different types of overcoating materials will allow different kinds of
20 testing to be carried out on samples.

21
22 In another embodiment of the present invention, light emitting diodes (LED) and
23 their related connection grids can be embedded into the slide surface itself. Using
24 surface mount LED sources or circuit mount laser sources it is possible to mount the
25 active light emitting source right in the slide itself and make the full system disposable
26 or reusable.

27
28 If the LED is mounted in the surface of the sample space and flush with it, it is
29 possible to put it right in the center of the image and then place the sample right on top
30 of the LED. If the sample is placed right on top of the LED the microscope objective

E0041280-032197

1 can be focused on the LED light bundle of rays and yielding an extremely bright
2 compact source right in the surface of the slide. If the LED spread function is matched
3 to the numerical aperture of the microscope objective it is possible to capture a very
4 large portion of the light from the LED and transfer it to the final image.

5
6 Because of the narrow spectral emission of LED's this method can be used to do
7 effective monochrome work, by choosing the appropriate LED emission range. For
8 instance using gallium nitride on silicon carbide LEDs it is possible to work in light
9 with a center wavelength 430 nanometers and which yields a very high resolution image
10 relatively free from chromatic aberrations.

11
12 In another embodiment of the present invention, a measuring photo diode can be
13 integrally mounted in the sample space to monitor the brightness of the illuminating
14 light and to act as a feed back system to stabilize the illuminating lamp supply. In this
15 way the photo diode is mounted very close to the target sample area and the photo diode
16 produces an electrical signal which is proportional to the illuminating light energy.
17 This electrical signal can then be used to control the power supply or to measure the
18 cumulative dose that the sample is being subject to.

19
20 Using different photo diodes, a clinician can monitor different spectral regions,
21 for instance in working with UV a silicon carbide or gallium nitride photo diode can
22 used. In monitoring visible light a silicon photo diode can be used, and in monitoring
23 infra red light an indium gallium arsenide photo diode can be used.

24
25 In another embodiment of the present invention, piezoelectric areas can be
26 provided in the slide, such as small microscopic piezoelectric transducers with their
27 related connection grids to create acoustic fields in the sample space to study the
28 reaction of acoustic waves on the sample. These same transducers can be used to detect
29 any acoustical emissions from the samples during examination and if these transducers
30 are used in the microperfusion version of these slides, described below, they can also

6041280-032197

- 1 provide information on flow rates and particulate loadings and fluid characteristics
2 using ultrasonic doppler techniques.
3

4 The packaging of a slide system in accordance with the present invention is
5 presently believed to be particularly advantageous. Conventionally, in the use of
6 microscope slides a microscopist has to clean the slide and then handle it carefully,
7 clean the cover glass and handle it carefully and then somehow combined the two
8 together with a suitable media and the sample without contaminating any of the internal
9 surfaces or the sample itself. This can be very difficult and is quite time consuming at
10 best. With the present invention, the slide will be prepared from a sterile package,
11 which is specifically made up with the slide and the cover glass located side by side in
12 the package and secured to the package sheet with a suitable removable adhesive
13 temporarily holding the slide and the cover glass in place and with a removable backing
14 sheet covering the adhesive area so that the adhesive does not stick to the rest of the
15 packaging material. The slide and the cover glass and the adhesive will all be clean,
16 contaminant free and sterile in their package. When the package is opened, the
17 microscopist can be assured that the slide is clean, the cover glass is clean and the
18 entire system was sterile.
19

20 The method of employing this system is to open the package in a manner similar
21 to the opening of a conventional bandage package, then remove the protective paper
22 from the adhesive and to then place the sample either on the cover glass or on the slide
23 itself and bring the two together and press them together tightly in order to activate the
24 adhesive.
25

26 The backing paper serves as a method of keeping the microscopist's hands away
27 from the cover glass and the slide while they are placing the sample in place and while
28 they are closing the slide. The backing paper also helps prevent the spread of any
29 hazardous material from the slide out onto the microscopist's hands. These slides in
30 their sterile packaging can either be supplied in single unit packages with one slide and

2004-03-20 09:13:37

1 one cover slip or they can be supplied in multiple unit packages in tape and reel form
2 n large reels with locating sprocket holes on the side to work in an automated slide
3 preparation system.

4
5 It is very important that no packaging adhesive remain on the bottom of the slide
6 within the sample area to disrupt microscope's light source. Accordingly, an adhesive
7 should be selected which will not leave an undesired residue and/or the packaging
8 adhesive may be located on the packaging materials such that it only contacts a
9 peripheral area of the slide. The sterility of the slide and the cover glass can be
10 accomplished either by radioactive sterilization or by heat sterilization prior to the
11 packaging. The benefit to the radioactive sterilization is it can be accomplished after
12 the slide has been packaged but, it is contemplated that in general the slide will
13 preferably have been packaged after thorough cleaning and in a clean room
14 atmosphere.

15
16 Depending on the nature of the backing material, and depending on its weight,
17 the backing paper can also serve as an element in an alignment system for automated
18 closing of slides or in a manual one-time package system by the provision of alignment
19 holes or sprocket holes. In particular, the backing material can include a score or other
20 indicia where the material will be folded to achieve alignment of the slide cover with
21 the sample area on the slide.

22
23 In another embodiment of the invention, the slide system can be employed to
24 create slides which function as micro perfusion chambers. In this embodiment, ports
25 can be provided either through the side of the slide or through the base of the slide into
26 the sample space to allow gases or nutrients or fluids or renewed samples to be
27 delivered to the sample space for examination or for assisting in the growth of an
28 organic or inorganic system or living system. In this way it is possible to watch the
29 formation of crystals in real time for instance in an inorganic system or it is possible to
30 watch a cellular, parasitic, bacterial or viral culture and keep it supplied with the

20041280-032197

1 nutrients that it requires over a long period of time. The flow through a microperfusion
2 chamber can also be used to provide cooling of the chamber.

3
4 It is also possible to use these micro perfusion chamber ports to supply different
5 substances which are desired to be tested. For example, in the area of cellular
6 research, it is possible to use these ports to supply different nutrients or different
7 chemical substances to observe the effect of chemicals on the cells. In this way it is
8 possible to use these slides to access the damage caused by chemicals which are desired
9 to test for biological hazards.

10

11 In this embodiment, the type of port can range from a simple drilled hole
12 through the base of the slide to an etched groove in the slide passing from the sample
13 chamber under the adhesive and out into the outside area. As an alternative, the ports
14 can be etched in a double layer base slide where two etchings are placed adjacent as the
15 slide is laminated together. The ports can be of several different styles and the ports
16 may number anywhere from a single port up to a large number of ports of different
17 sizes depending on their desired use.

18

19 As an example, in a complex microperfusion chamber there can be a pair of
20 main ports for incoming and outgoing nutrient, a set of smaller ports for incoming and
21 outgoing fluids, and a set of even smaller ports for incoming and outgoing gases. On
22 the microscope base that the slide would fit into, there will be a port sealing adapter
23 which will involve ports and suitable sealing devices, such as O rings or lip seals, to
24 seal the ports from leaking material or gases from the slide into the surroundings.
25 These ports can either be disposable and made of a cast material or they could be
26 permanent, cleanable and sterilizable ports which could be removed from the slide
27 holding base. Connections to the disposable or reusable port adapter can be made via
28 hoses to the various supplies of nutrients, chemicals, samples or gases and their related
29 recirculating devices.

30

60041260.0324099

1 It is important in using these slides to deliver a precisely measured quantity of
2 material to the sample space to insure that one doesn't overflow the expansion wells or
3 put so much material that the function of the expansion wells would be compromised
4 resulting in the adhesive rings being breached. This can be accomplished using
5 pipetting or micro pipetting techniques with disposable pipettes. The pipettes will be
6 calibrated to work with the slide system and deliver a precisely measured quantity of
7 sample to the sample space so that there was little hazard of the sample being excessive
8 and causing damage to the system or leaking from the system.

9
10 As will be apparent to those of skill in the art, the present invention provides a
11 range of advantages and improvements over conventional slides for microscopy.
12 Embodiments of the present invention can be used for live or fixed specimens of
13 biological, organic or inorganic matter. The present invention is believed to be suited
14 for use with living samples of cells, vital fluids such as blood and lymphatic fluid,
15 parasites, bacteria or virus or combinations of these. The present invention can also be
16 used to study the interaction of living material with chemicals; biological products;
17 electric, magnetic, photo, acoustic or ionizing radiations; and other living material.

18
19 Brief Description of the Drawings

20
21 The most basic form of the invention is illustrated in figure 1. A rectangular
22 slide of standard form (A) has a ring of adhesive (H) coated onto one surface. The ring
23 of adhesive may be controlled in thickness to provide the desired sample thickness. An
24 amount of sample material is placed in the center of the slide at location (D). A
25 coverglass (C) is applied to the adhesive ring and pressed into place to spread the
26 sample material over a portion of the area contained by the adhesive ring and to activate
27 the adhesive.

28
29 In Figure 2 the same basic slide included an expansion zone or moat shown as
30 (E) inside the adhesive ring (D).

50041290.03P197

1
2 Figure 3 shows the addition of an antiseptic or antibiotic ring (F) inside the
3 adhesive ring (B) and outside the expansion zone (E).

4
5 Any of the slide systems described herein can also be supplied in clean and
6 sterile form by providing both the coverglass and the previously adhesive coated slide
7 in a sealed paper or plastic film protective package. In Figure 4 the slide (A) and the
8 coverglass (C) are contained in a paper wrapper. The lower layer (G) of the paper
9 wrapper holds the slide and the coverglass in place with areas of low tack adhesive
10 which peel away from the slide and coverglass freely leaving no optical residue once
11 the slide is prepared. An upper layer of the paper has a release layer coated onto the
12 inner surface so that it will peel away from the adhesive ring (B) without causing any
13 loss of adhesion. To use the system, the upper paper is peeled off the lower paper (G)
14 to reveal the slide and the coverglass. The sample material is applied to the slide or the
15 coverglass in locations (D) or (V) respectively. The lower paper is folded along the
16 marked or scored fold line (I) which ensures that the cover glass (C) will be positioned
17 properly on the adhesive ring (B). Once the coverglass has contacted the adhesive ring
18 pressure is applied to the sandwich to set the adhesive. The paper (G), which is now on
19 the top and bottom of the sandwich of layers, is peeled away to leave the finished slide.
20

21 In automated slide preparation systems it is desirable to provide the packaged
22 slides and coverglasses in continuous tape form with locating holes for alignment in an
23 assembling machine. The tape can be supplied in reel or zigzag stored form.
24

25 In Figure 5, the lower paper is provided with sprocket holes (J) on both outer
26 boundaries to allow the slides to be aligned during assembly. The upper paper strip is
27 removed from the lower tape and stored on a roller, the lower tape is brought into the
28 area of a sample dispenser which delivers a measured amount of the sample material to
29 the location (D) or (C) of the slide. The lower paper (G) is then folded along fold line
30 (I) by the carrier to provide the preliminary seal. The lower paper is then peeled away

60041230-032197

1 from the finished slide and the slide is delivered to a final closing press which closes
2 the coverglass to the slide with the required final closing tolerance. The closing height
3 is then checked and when confirmed the slide is placed in a delivery tray. If the slide is
4 not closed to the required tolerance it is rejected and placed in a reject tray. The slide
5 is labeled with a numerical label to identify the slide as to contents.

6
7 Figure 6 shows a version of the slide where dots of a test material (K) have been
8 coated or printed onto the base slide (A). This test material could be any material
9 which might interact with the sample such as an antiseptic, antibiotic, chemical,
10 antibody, stain, fluorochrome or other living organism.

11
12 Figure 7 shows a version of the slide where expansion zones are formed in the
13 base slide (A) in the form of wells recessed into the slide surface.

14
15 Figure 8 shows a slide for testing of electrical properties of a samples such as
16 electrophoresis, electrolysis, electrochemical action, corrosion or dipolar behaviour.
17 Conductive traces (R) are coated onto or etched onto the base slide (A) using
18 semiconductor fabrication or printed circuit board techniques. These conductors are
19 connected to terminals outside the sample space (Q) for connection to electrical sources.
20 Conductive circles (S) prevent field concentrations at the ends of the conductors.

21
22 Figure 9 shows the same concept as Figure 8 except the conductors are coated
23 with a material (T) which may be a dielectric, an insulator, a semiconductor, a
24 reflective layer or layers, or other material.

25
26 Figure 10 shows a version of the slide shown in Figure 8 in which the conductor
27 (R) and (R') continue through the sample space and terminate at a connection terminal
28 (Q') and (Q''). This allows current to flow along the conductor creating a magnetic
29 field in the sample chamber. The use of conductors (R) and (R') means that both

60041280.032197

1 magnetic fields and electrical fields can be established in the sample space at the same
2 time.

3

4 Figure 11 shows a version of Figure 10 in which the conductors are isolated
5 from the sample space by a layer of material (U) as in Figure 9 described above.

6

7 Figure 12 shows a slide in which the base slide has been marked with a
8 reference scale or reticule (O) and with a serial number (P).

50041280.032137

1 I claim:

2

3 1. A slide system for microscopy comprising:

4 a slide including a sample holding area enclosed by an adhesive area; and

5 a slide cover to engage said adhesive area such that said sample holding area is
6 sealingly covered by said slide cover and defines a desired resulting total thickness of
7 said sample holding area.

8

9 2. A method of preparing a slide for microscopy, comprising the steps of:

10 (i) placing a sample on a slide within a sample holding area enclosed by an
11 adhesive material;

12 (ii) locating a slide cover over said sample holding area to engage said adhesive
13 material; and

14 (iii) pressing said slide cover to said adhesive material such that said sample
15 holding area is sealed by said adhesive and to obtain a desired thickness for said sample
16 area.

17

267280 08214005

Abstract

1
2
3
4
5 A system of components designed to allow the quick, sterile and safe
6 preparation, use and destruction of microscope slides containing living or hazardous
7 preparations. The system consists of a slide and coverglass which are precleaned and
8 supplied in a sterile wrapper designed to facilitate easy handling. The coverglass or the
9 slide have an outer ring of adhesive coating which may include a built in antiseptic
10 compound. In one embodiment the slide has one or two moats engraved in the surface,
11 the first moat to accept excess material and the second moat to isolate the sample from
12 the antiseptic adhesive. The slide and coverglass, when closed have substantially exact
13 dimensional tolerances to allow them to be placed in holders on a microscope stage.
14 When used in appropriate holders and engaged to the microscope objective, the objective
15 will be focussed in the center of the z axis of the sample space. The slides can include
16 a number of different additional features to allow live cell and culture work and to
17 provide varying environments in the slide sample chamber.

60041260:032197

60/041280

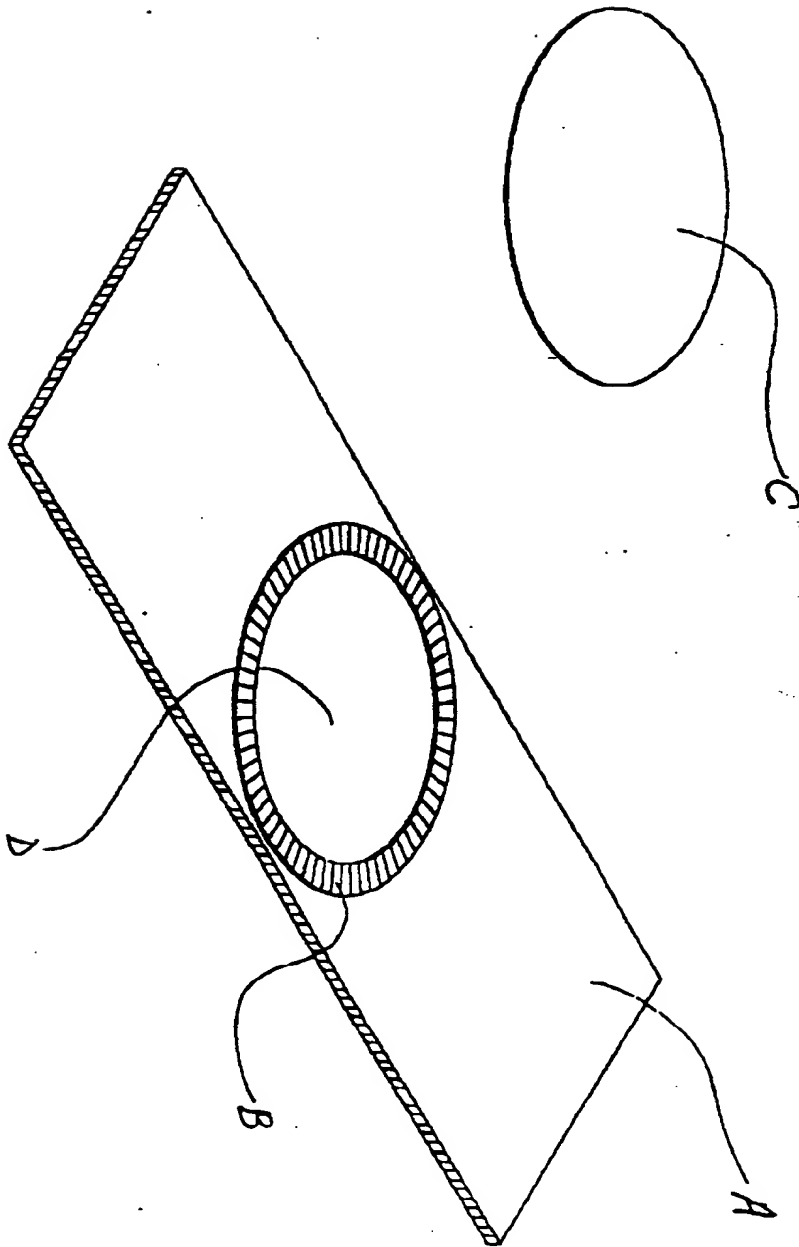


FIGURE 1

60041280.032197

60/041280

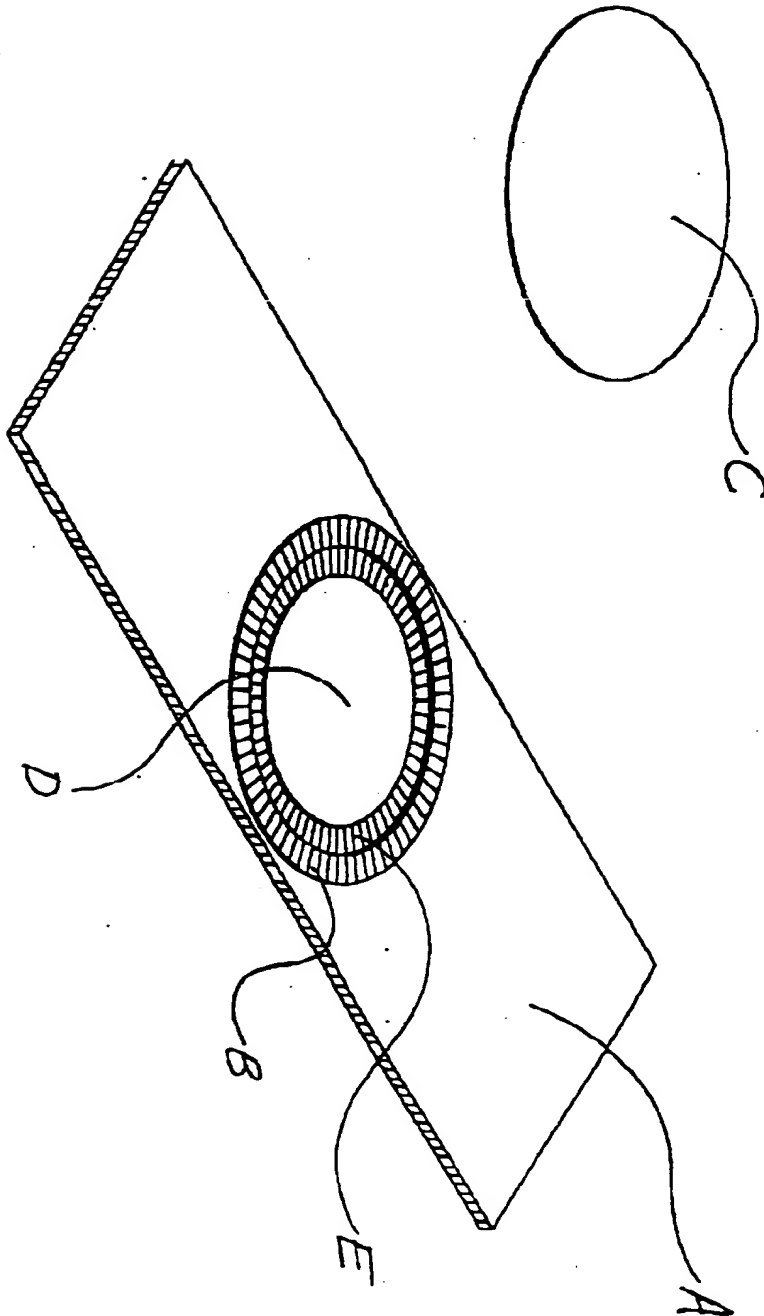


FIGURE 2

60041280.026197

60/041280

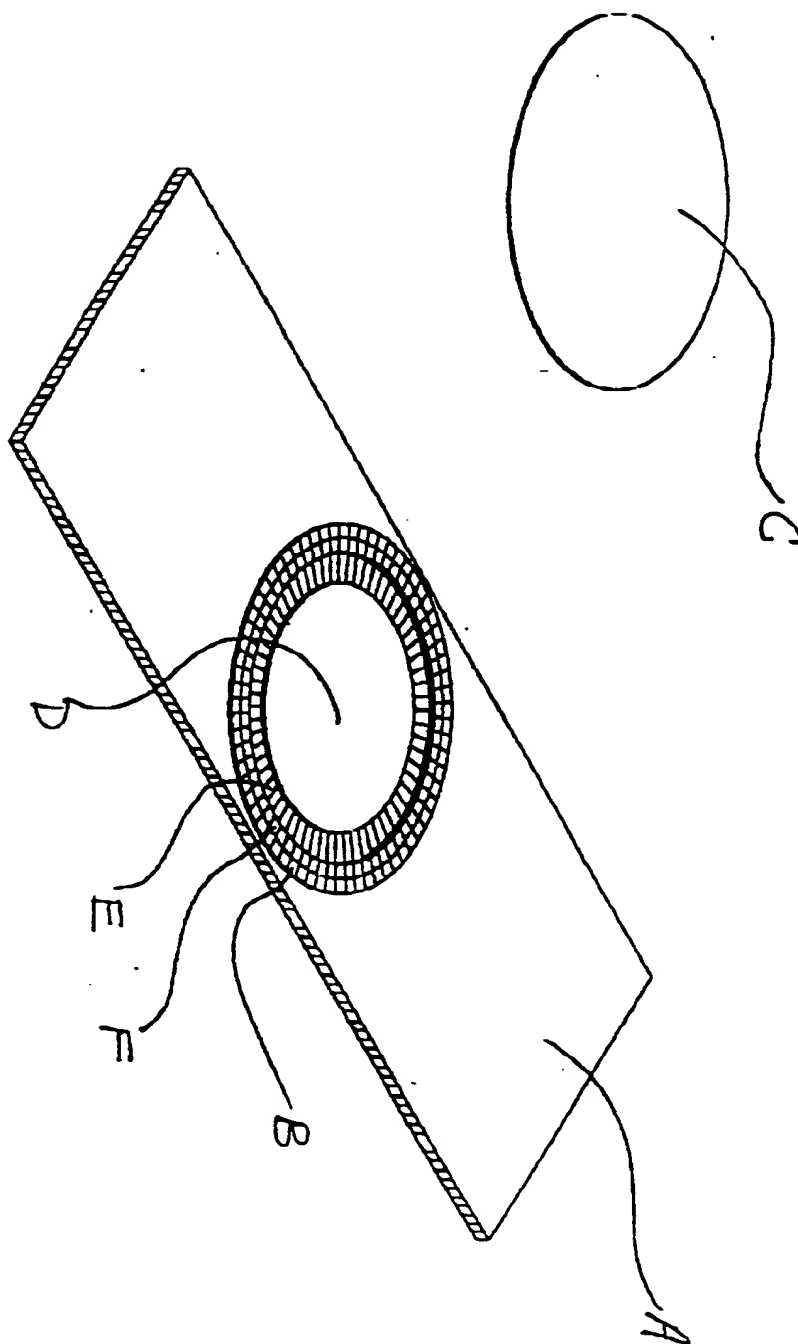


FIGURE 3

60041280,032197

60/041280

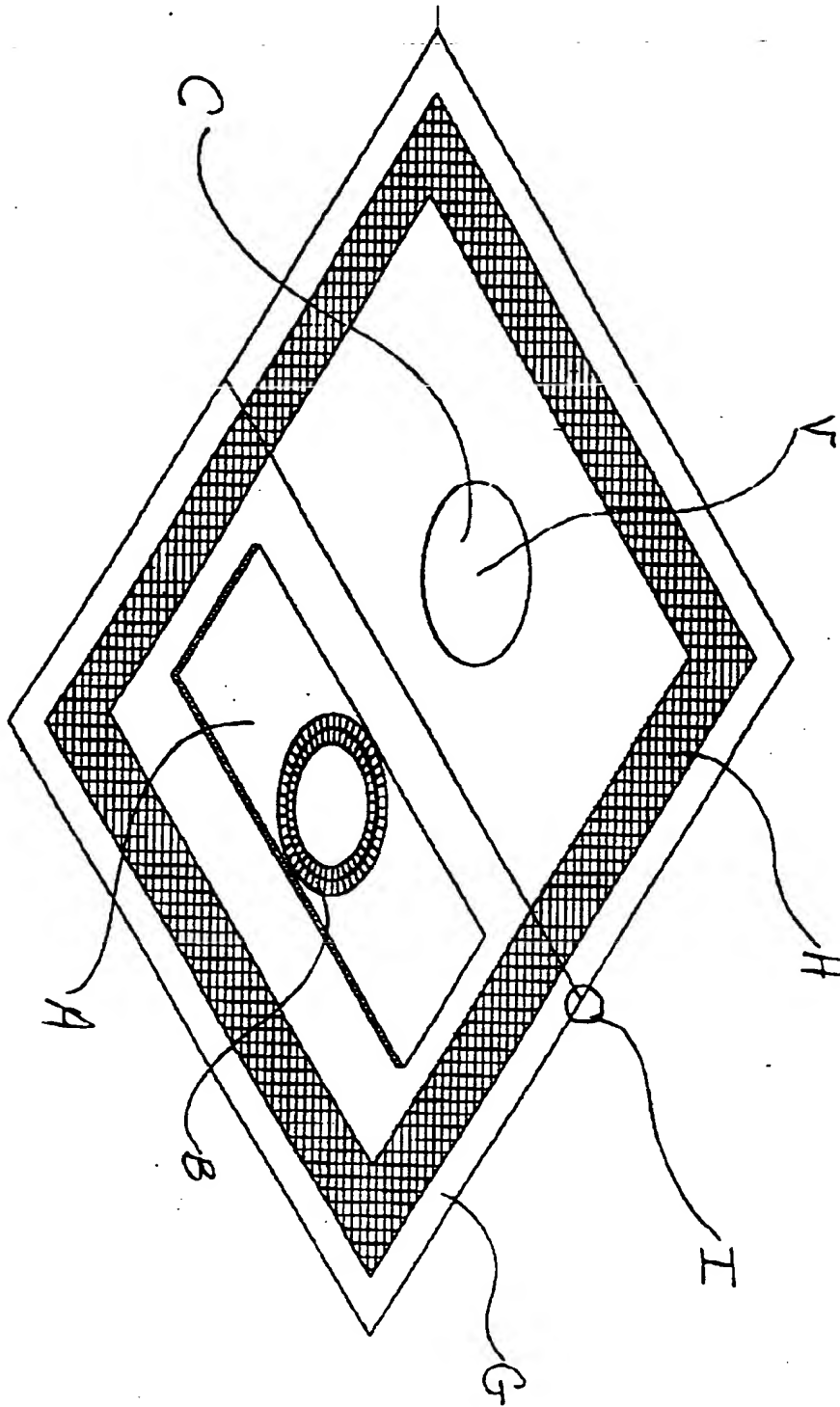


FIGURE 4

50041280.032197

60/041280

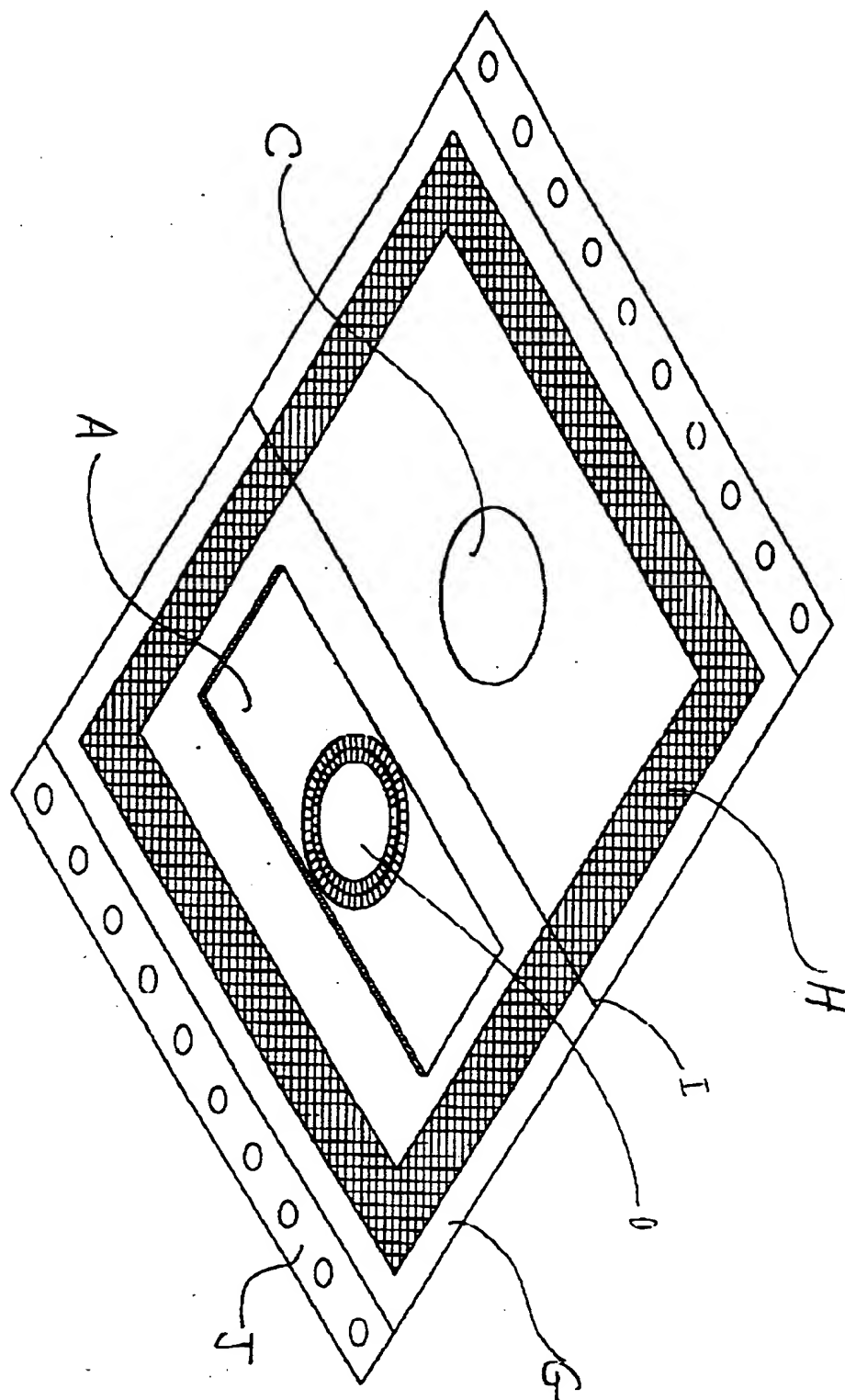


FIGURE 5

60041280.032197

60/041280

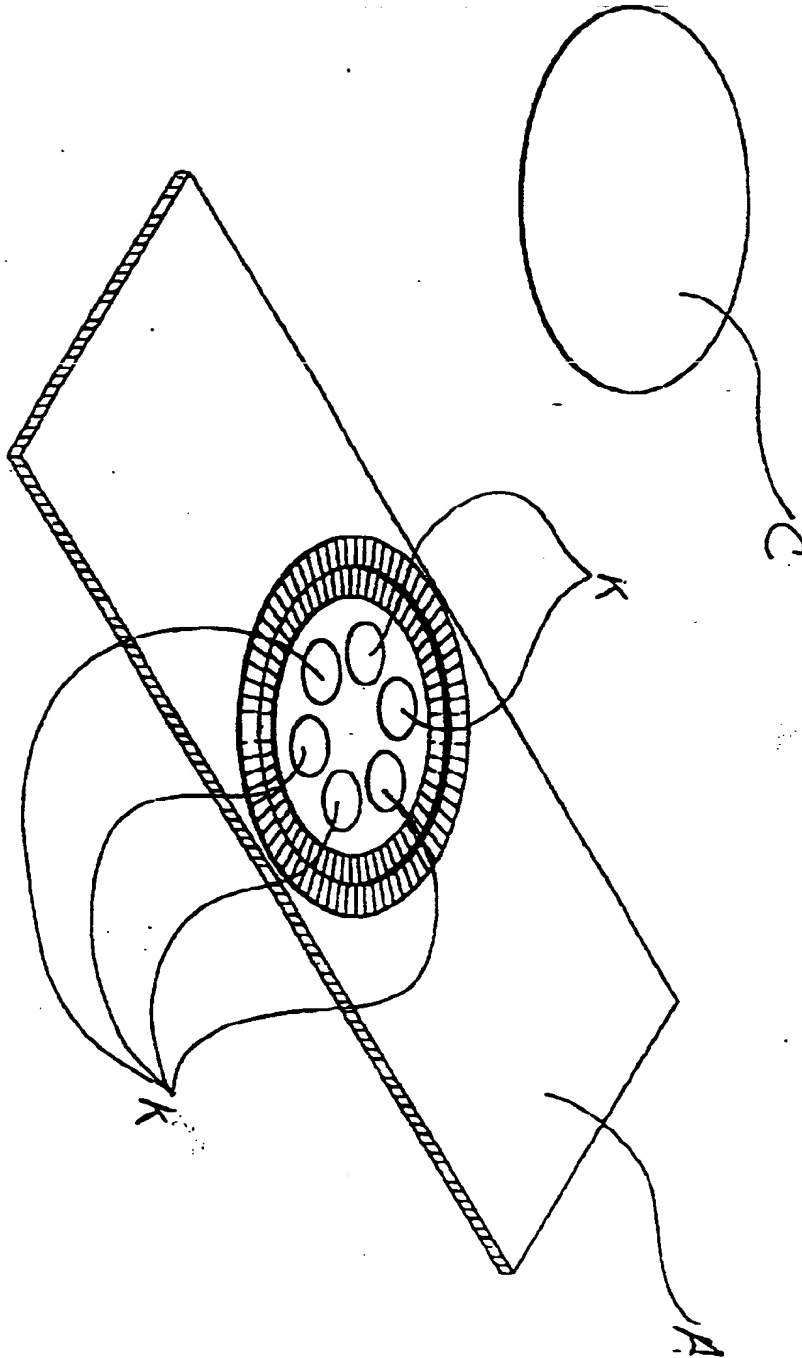


FIGURE 6

60041280.032197

60/041280

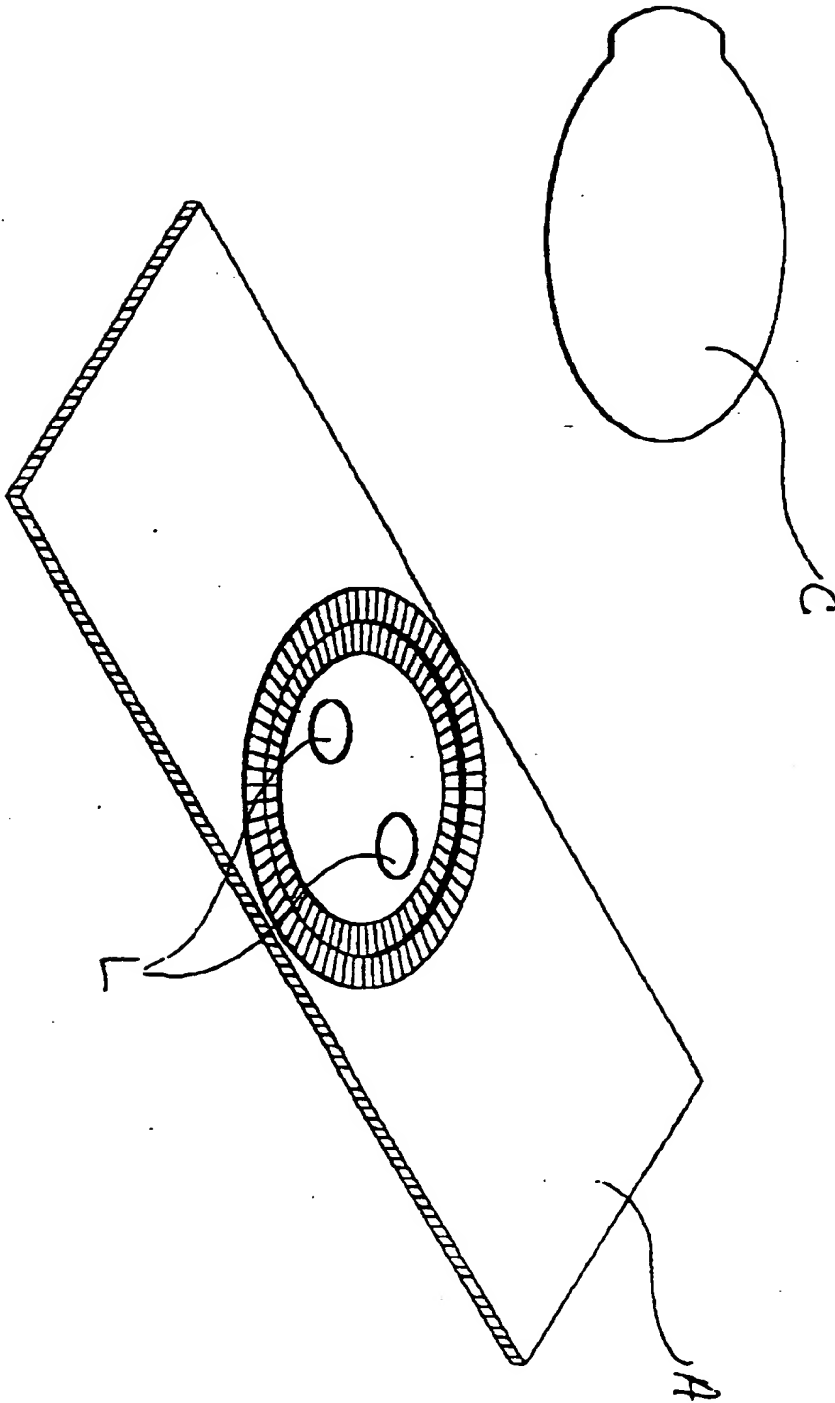


FIGURE 7

60041280.032197

60/041280

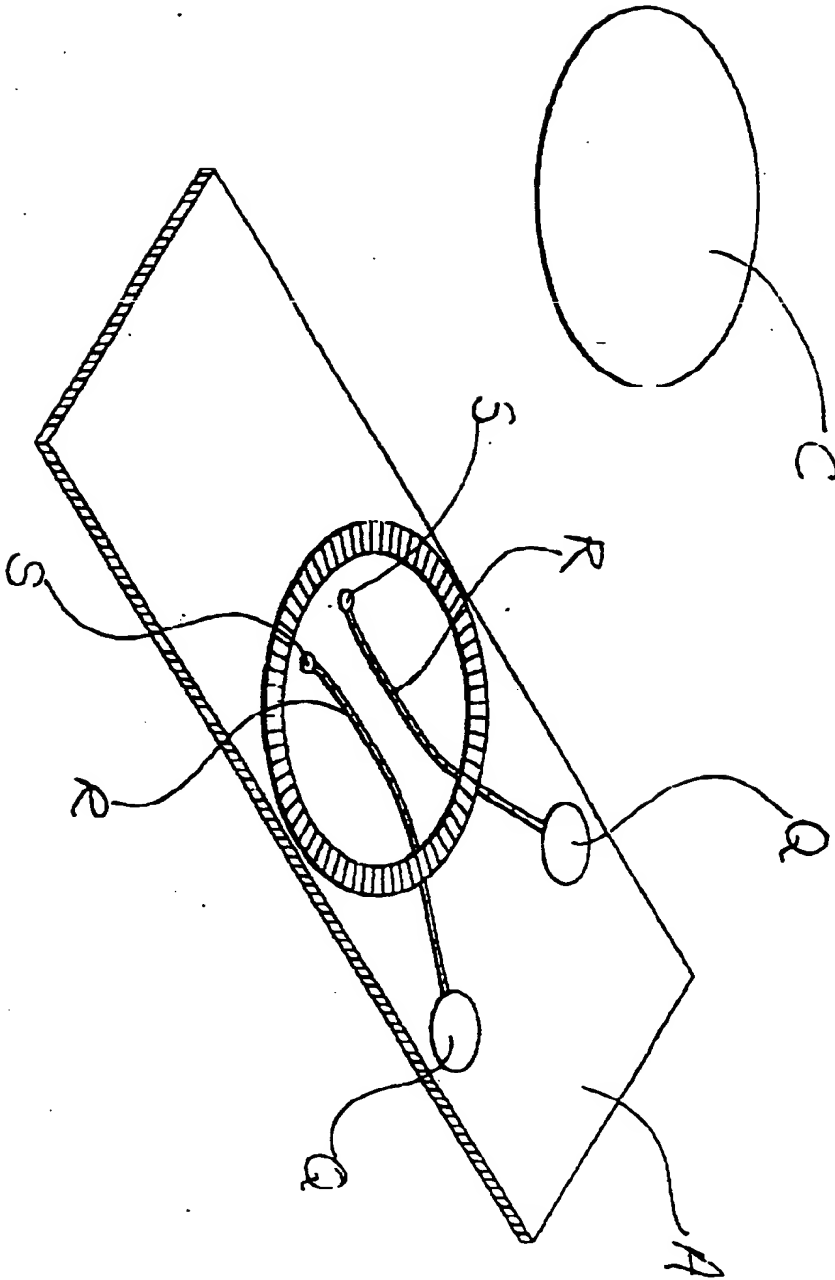


FIGURE 8

60041280-022197

FIGURE 8

60/041280

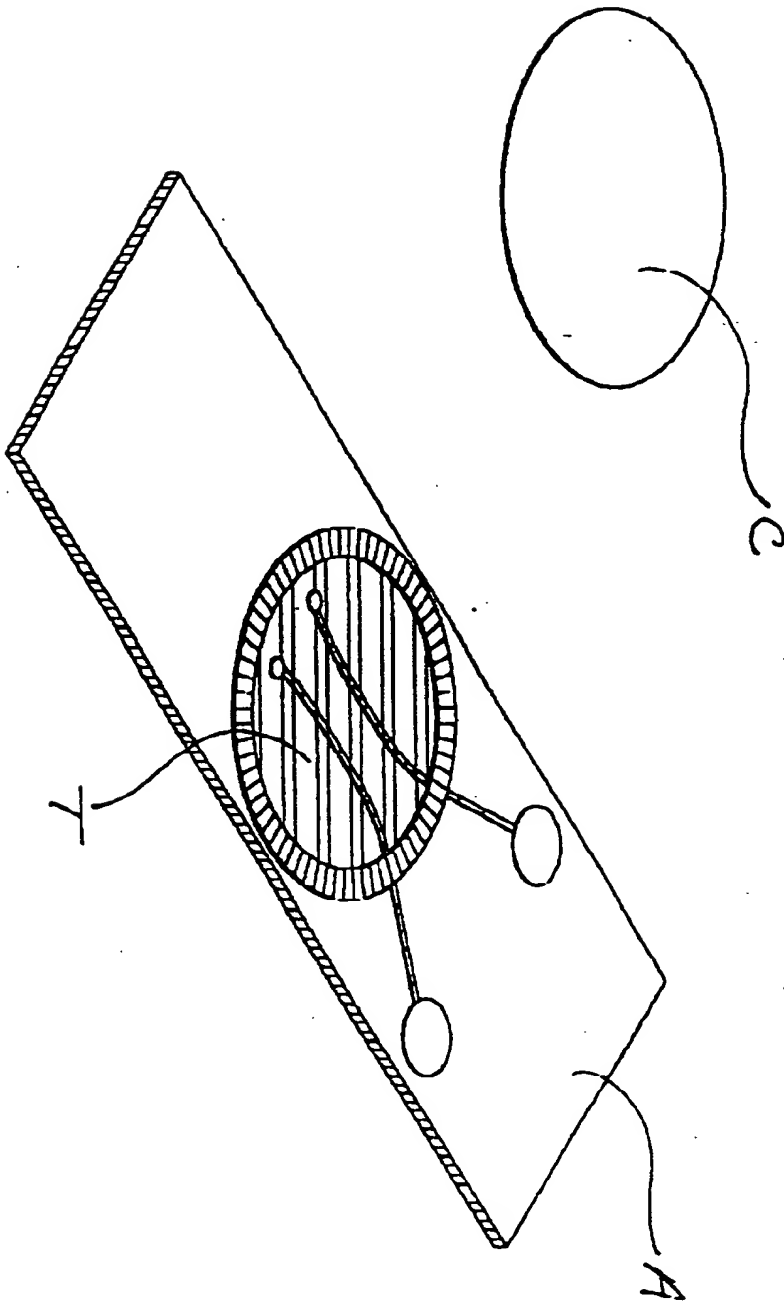


FIGURE 9

60041280.032197

60/041280

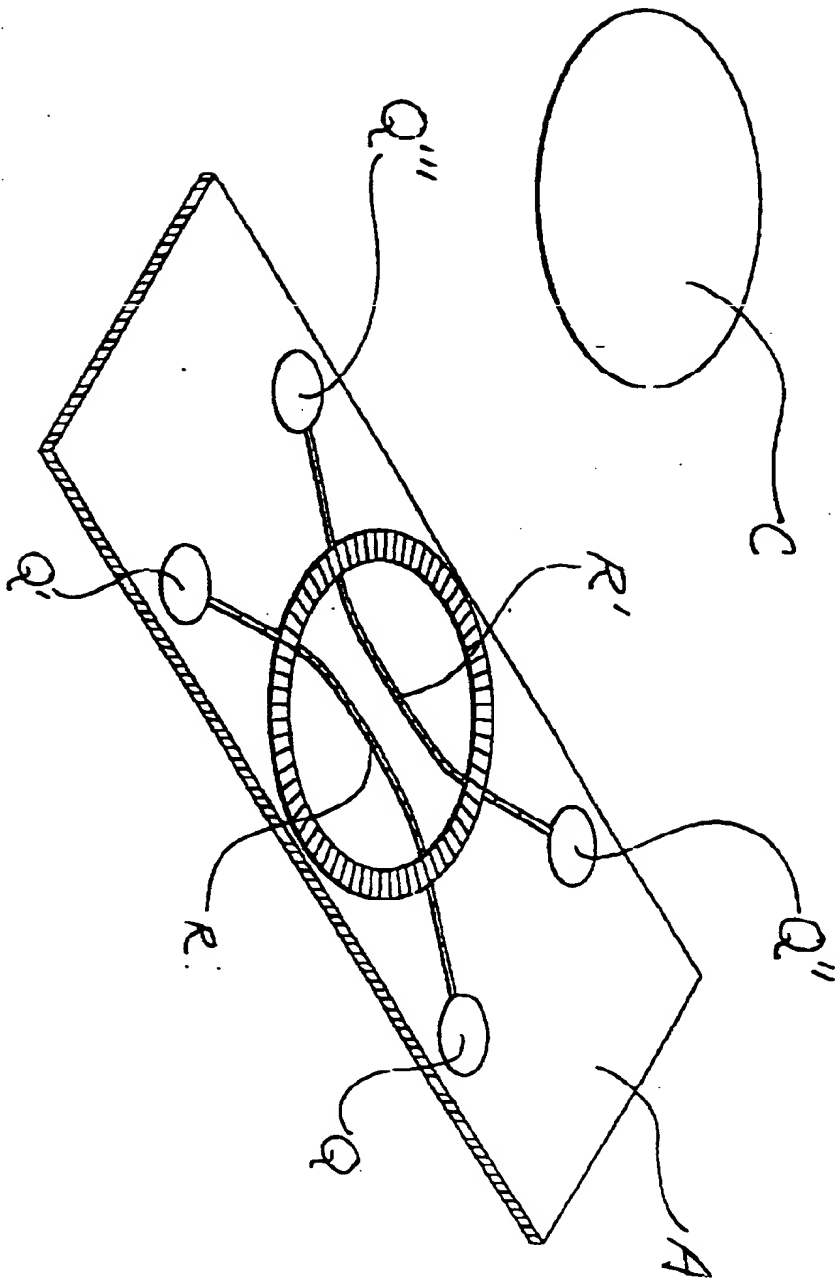


FIGURE 10

60041280.032197

60/041280

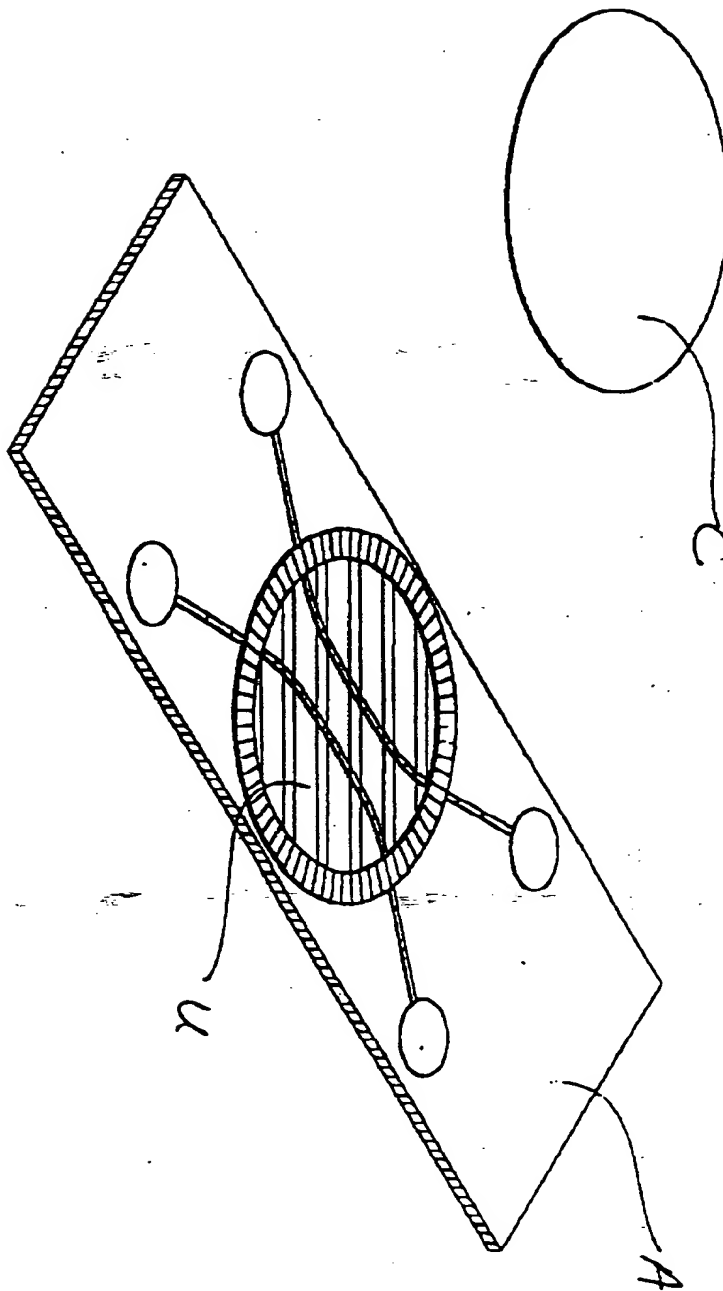


FIGURE 11

60041280.032157

60/041280

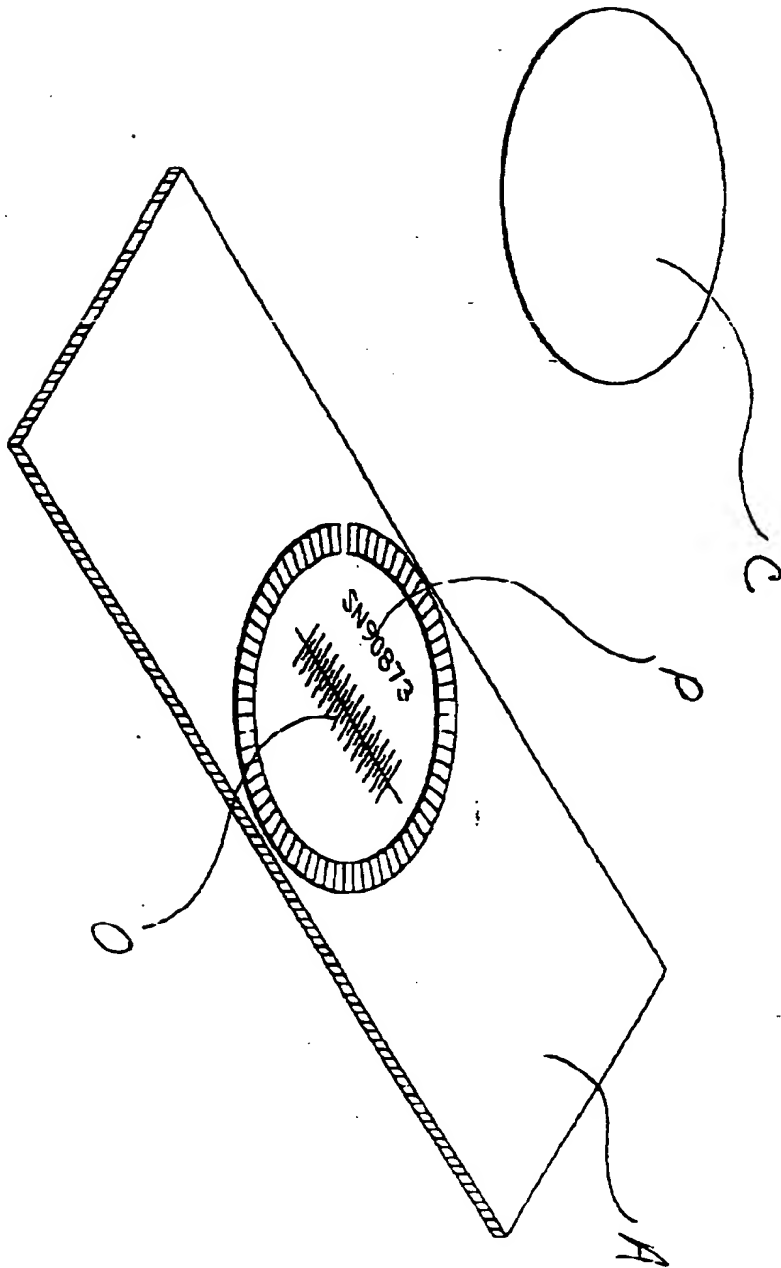


FIGURE 12

60041280.03E137